



# Jordco® e-Foam® Guide

## Endo Foam and High Density Inserts


Jordco's e-Foam insert is the standard for all clinicians who perform endodontic procedures. Our patented foam inserts are designed to be used with all Jordco EndoRings® and they provide optimum performance for superior cleaning, storing and organizing of endodontic files.

### e-Foam® Rotary HD

High Density Foam Inserts

NEW COLOR  
Jordco Silver

ERFSHD-s

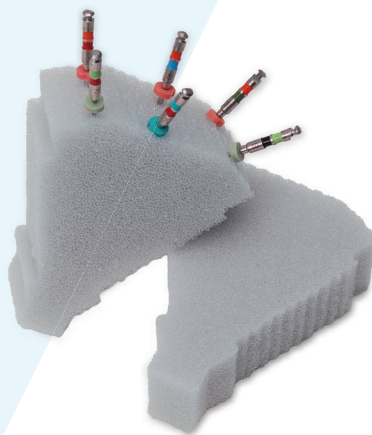


Jordco e-Foam Rotary HD (High Density) inserts are specially designed for cleaning rotary instruments. A rotating file secured in a hand piece can be placed directly into the high density foam to assist in cleaning debris from rotary flutes. Our Rotary HD foam minimizes tearing or shearing of the foam while a file is in motion.

### e-Foam®

Regular Density  
Foam Inserts

TWO COLORS

ERF-s (Jordco Blue & Jordco Yellow)  
ERFB-s (Jordco Blue)  
ERFY-s (Jordco Yellow)

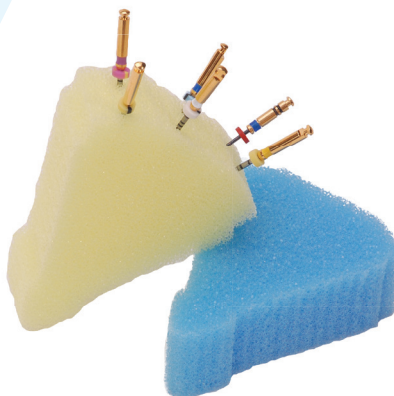

Single use patented open cell yellow and blue foam. Can be steam autoclave processed prior to use. During steam autoclave processing, endodontic instruments can be inserted into the foam. Jordco's e-Foam is available in 2 colors to allow the option of color coding instrument setups.



## FEATURES AND BENEFITS | ALL FOAM INSERTS

### Superior Design and Cleaning

The e-Foam inserts cell structure has a much higher and desirable firmness than most all of the imported imitators. Jordco patented e-Foam efficiently stores and cleans debris from endodontic files. All of our foam is designed and made in the USA.

### Three Colors of e-Foam Inserts

Jordco blue and yellow e-Foam inserts, and silver HD Foam inserts, provide additional options for color coding instrument set-ups.

### Locking Tail Feature

The shape of the tail portion of all the foam inserts has a positive locking feature to engage the EndoRing cup and prevent dislodgement while in use. The unique shape of the tail also helps to identify the Jordco brand e-Foam insert.

### Foam Can be Autoclave Processed

All Jordco foam inserts can be loaded with new or clean files and steam autoclave processed before a procedure. The open cells structure allows for easy penetration of steam.

### Fits All Jordco EndoRings

All Jordco foam inserts are designed to fit both the Original EndoRing® and the EndoRing® II products.

### Single Use

Single use foam inserts help prevent cross contamination and prevent accidental sticks. Used files can be placed into the foam insert for disposal in a sharps disposal container.

### Faceted Working Surface

The top of the foam insert has three distinct facets that help organize files for ease of use during a procedure. Different tapers or diameter of files can be easily identified and located during a procedure.

## UNIQUE FEATURES AND BENEFITS | E-FOAM ROTARY HD

### Perfect for Rotary Files

Jordco's Rotary HD (High Density) foam is specially designed with a higher shear strength than our regular e-Foam. This helps to minimize the tearing of the foam while the rotary file is in motion and being inserted into the foam.

### Superior Cleaning

While in motion, rotary files can be inserted into the HD Foam to rapidly and efficiently remove both dentin and pulpal debris.