




# EndoRing® II

## Hand-held Endodontic Instrument With Reduced Glare Metal Ruler

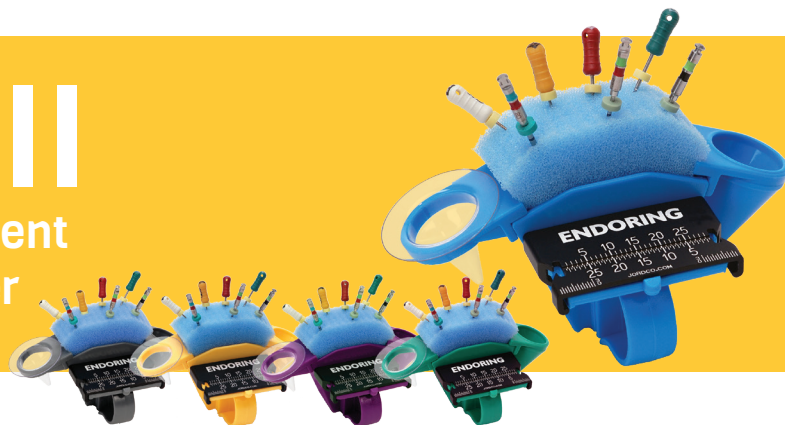

### EndoRing II Endodontic Instrument Family of Products

The EndoRing II organizer enables the clinician to directly place, store, measure and clean endodontic hand and rotary instruments within the operating field. It reduces the risk of cross-contamination and minimizes the passing of sharps between dentists and staff. Using the EndoRing II organizer helps make endodontic procedures safer and faster.

### TECHNICAL SUMMARY

**EndoRing II:** Assemblies are molded from an autoclavable blue thermoplastic resin. Assembly consists of six pieces: 1 EndoRing Cup, 3 Finger Rings (S, M, L), 1 reusable GelWell Cup, and 1 GelWell Adaptor Ring. Note: the EndoRing Cup has a molded plastic 30-mm scale for right or left handed users.

**EndoRing Metal Ruler:** Autoclavable anodized black aluminum detachable ruler with contrasting white laser-etched scales in 1mm and 1/2mm increments. Accessory scale provides marking up to 37mm for extra long canals. Ruler can be used by right and left handed users. Available with rubber stop locks.

**GelWell System:** Reusable Cup: 2ml cup is molded from an autoclavable blue thermoplastic resin and attaches to the EndoRing II organizer. Single Use Cups: 0.65ml cups molded from a clear thermoplastic resin. The single use cup is placed into the holder and attaches to the EndoRing II organizer.

**e-Foam:** Single use patented open cell yellow and blue foam. Can be steam autoclave processed prior to use. During steam autoclave processing, endodontic instruments can be inserted into the foam.

### PRODUCT SPECS

**Name:** EndoRing II with Metal Ruler  
**REF:** ER2RG-s (Blue), ER2GRG-s (Green),  
ER2YRG-s (Yellow), ER2VRG-s (Violet),  
ER2GYRG-s (Gray)  
**Weight:** 5.5 oz

### PACKAGE CONTENTS

- 1 EndoRing II Assembly
- 4 Single Use GelWell Cups
- 6 e-Foam Inserts
- 1 Black Metal Non Glare Ruler

### Each Assembly Includes:

- 1 EndoRing Cup
- 3 Finger Rings (S, M, L)
- 1 Reusable GelWell Cup
- 1 GelWell Adaptor Ring

### FEATURES AND BENEFITS

#### Hand-held Instrument

Allows direct operator control of hand and rotary instruments. Eliminates time consuming passing of dangerous files from assistant or bracket tray. Minimizes eyestrain and loss of attention from searching for files, etc. Saves several minutes each procedure.

#### Held on Index Finger

Frees both hands. Allows operator to hold mouth mirror in one hand while selecting and operating files or instruments with the opposite hand. Includes three ring sizes (small, medium and large). The small ring is ideal for the petite hand.

#### Two Color EndoFoam Cleans and Organizes

Jordco blue and yellow e-Foam inserts provide an additional option for color coding instrument set-ups. While in use, e-Foam efficiently cleans bioburden from most all endodontic files.

#### Detachable Laser Etched Metal Ruler

Autoclavable metal ruler with a precise laser etched right or left-handed 30mm scale. Includes 1/2 mm calibration markings on the lower half of the scale. The metal ruler also includes a secondary scale to measure files up to 37-mm in length. Includes a calibration track for easy positioning of instruments. Ruler can be easily detached for cleaning and processing.

#### Easy Calibration of Curved Instruments

Curved files can be easily positioned into the calibrating track within the standard or metal ruler and temporarily straightened with slight pressure. After the stop is adjusted, the file is removed and allowed to reassume its curved shape.



Jordco, Incorporated  
8531 NE Cornell Rd, STE 800, Hillsboro, Oregon 97124 USA  
(800) 752-2812 toll free USA (503) 531-3904 www.jordco.com



U.S. Patents D689190, 6036490, 8231734, D659837, 8231734B2, 8635735B2

REV2209